

Recent Technical Developments in the Standardized Separation and Measurement of Extracellular Vesicles

Amy Phillips¹, Robert Vogel², and Murray F. Broom³

¹Izon Science US Ltd, Christchurch, New Zealand (amy.phillips@izon.com), ²The University of Queensland, Brisbane, QLD, Australia, ³Izon Science Ltd, Christchurch, New Zealand

www.izon.com

The Need for Standardised Analysis In Extracellular Vesicle Research

The concentration, size, and charge of vesicles reflect the health of the organism or tissue they are extracted from¹. Accurate measurement of these parameters, therefore, has critical implications in the clinic. Furthermore, sample preparation can influence the vesicle characteristics².

A standardised method of vesicle preparation and analysis is needed to allow comparison between research groups.


Current purification methods have problems with aggregation, or long, laborious processes. Commonly used characterization techniques lack sensitivity, and resolution³.

We have developed a method of isolation and characterization, which together provide a sensitive method for determining vesicle size and concentration.



TRPS: Detailed Size, Charge and Concentration of EVs

Izon TRPS use a tunable elastic pore to measure particle properties. The size, charge and concentration of vesicles are determined with high resolution and sensitivity from the ionic current pulse signal they generate as they pass through the elastic pore⁴.

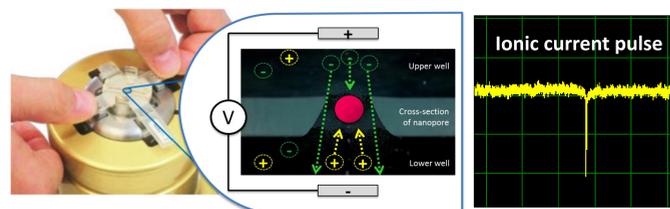

Tuning the pore size, applied pressure and voltage improves the measurement sensitivity and enables a wider range of particles to be analysed.


Particle Size: 50nm – 20µm+
Conc. Range: 10⁵ – 10¹² particles/mL
Sample Size: 30µL
Sample Media: Aqueous Buffer

Size-Based Concentration

A new definition of particle concentration based on the size range being examined has been proposed, e.g. C₇₀₋₂₀₀



This is more meaningful than the standard definition in that it specifies the range under consideration and of interest to the researchers (e.g. 70 – 200nm), and will aid comparability between measurements.

TRPS is unique in its ability to simultaneously determine the size and concentration (particles/mL) of vesicle samples based on the rate at which they pass through the pore.

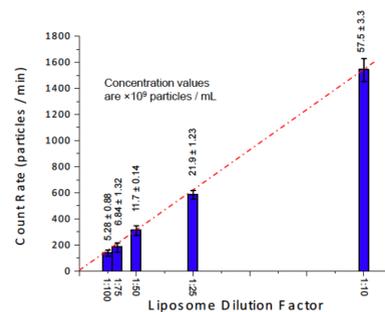


EV Identification

Aptamer based tool kit in development, to be used in conjunction with qNano to provide biological information.

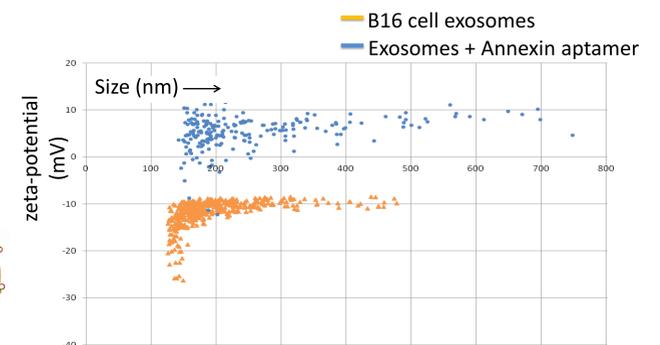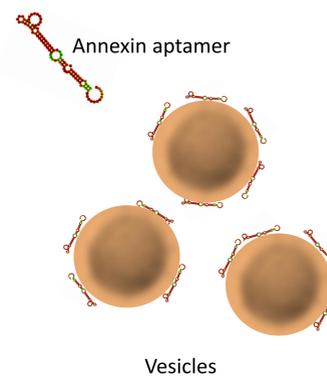


Particle Surface

Analysis of ζ-potential analysis & surface modification

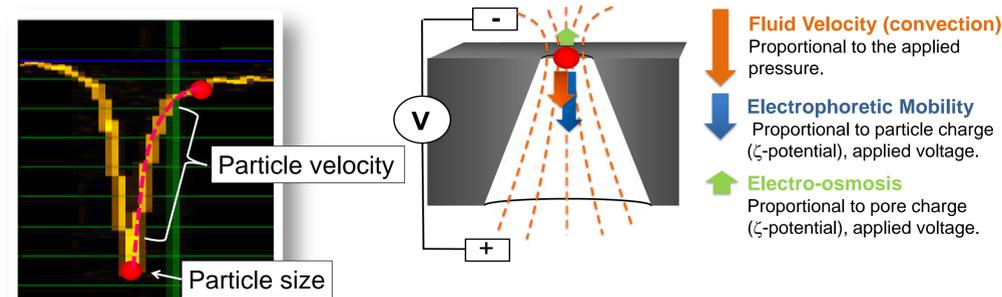

ζ-potential is calculated from the shape of the resistive pulse signal generated. Changes to the particle surface through, for example, protein or aptamer binding can cause changes in ζ-potential. Due to the particle-by-particle measurement, mixtures of vesicles with different surface properties can be measured.


Surface charge of vesicles decreases vesicles protein binding. Ratios of differently charge vesicles can be detected.

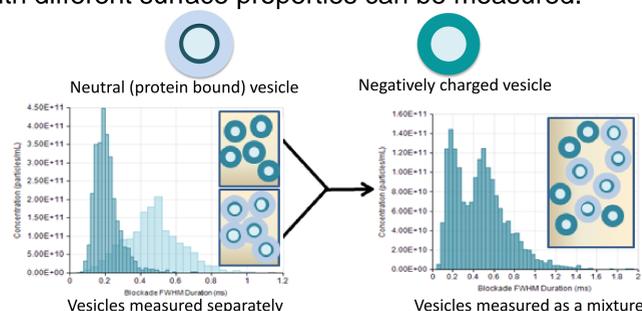

The interaction of aptamers with biological targets on the vesicle surface, e.g. CD63, can provide discriminatory information based on the size and charge profiles of the different populations.

Conclusions

- TRPS is the only high resolution technique which measures Size and Count, as well as surface Biomarkers (Figure)
- Size-Based Concentration is essential for meaningful comparison of measurements
- qEV columns provide a simplified approach to EV separation and isolation, and are expected to become the standard method.

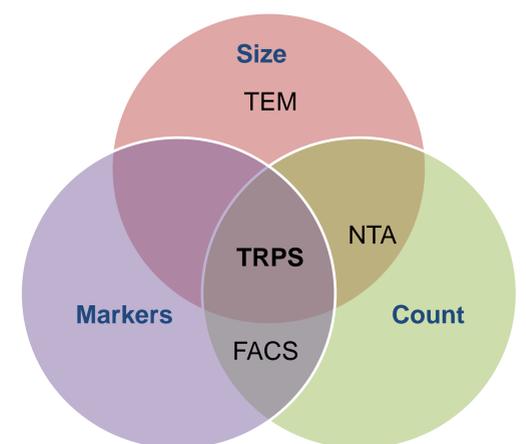

Comparison of the capability of techniques for analysing vesicles.

References

- van der Pol, *et al.*, (2012) Classification, functions, and clinical relevance of extracellular vesicles, *Pharmacol Rev* 64, 1-33
- Sowery *et al.*, (2007) Dynamically resizable nanometer-scale apertures for molecular sensing, *Sensors & Act B* 123, 325-330
- Yuana *et al.*, (2011) Pre-analytical and analytical issues in the analysis of blood microparticles, *Thromb Haemost* 105, 396-408
- Anderson *et al.*, (2013) A comparative study of submicron particle sizing platforms: Accuracy, precision and resolution analysis of polydisperse particle size distributions, *J Coll Inter Sci* 405, 322-330
- Roberts *et al.*, (2012) Biosensors and Bioelectronics 31, 17-25